




Truly interoperable connectivity

Mindray's flexible connectivity options enable you to future-proof your hospitals systems, for scalable growth and future innovation. Our vendor neutral and comprehensive connectivity solutions make digital patient data integration simple and secure.

- Seamlessly connects to your EMR systems
- Enabled caregiver login and patient identification
- Hybrid BeneVision central monitoring system to maximize your investment
- M-IoT device manager provided to screen the clinical efficiency of the equipment

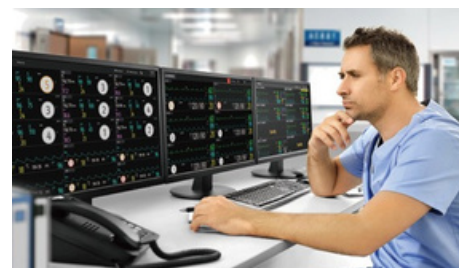

VS 9

Vital Signs Monitor

True to Life

True to form quality and reliability

Created with longevity in mind, the VS 9 is supported by innovative maintenance software and a robust design.

Trouble-free Cleaning

Validated with up to 51 leading cleaning agents, to aid faster wipe-downs and less degradation over time.

Easier installation

Step-by-step configured for your own clinical requirements with a smartphone-style Startup Wizard.



Useful workboard with optional extended battery for longer runtime

Tool-free design for faster installation under one minute

Dedicated accessory management as you prefer

Quiet and smooth casters for ease of use



Request a quote

Sales / customer service
sales@norsomedical.co.uk (NI)
sales@norsomedical.ie (ROI)

Website

www.norsomedical.com

Phone Number

Call 9:00 a.m. - 5:00 p.m. Monday - Friday
02890 343927 (NI)
01 5314337 (ROI)



See the true picture of patient deterioration

To help ensure a timely reaction to patient deterioration and greater compliance with hospital guidelines, the VS 9 provides a comprehensive range of vital signs parameters with a variety of tools to support clinical decisions.

Automatic EWS

- Fully compliant with all EWS standards
- Provide customizable escalation messages

Orthostatic BP

- Innovatively help to assess and avoid fall risk
- Customizable due to your own clinical requirements

BP Average

- Practice unattended AOBP by your own workflow
- Provide accurate values for real patient status

From inpatient admissions to spot checks, the VS 9 Vital Signs Monitor unlocks the full potential of routine observations with comprehensive and flexible workflow options empowering clinicians to deliver high quality bedside care to each and every patient.

True-to-life vital signs

The VS 9 delivers precision vital signs measurements at rapid speeds, maximizing patient comfort and clinical efficiency.

Mindray TrueBPTM inflation NIBP measurement in 15 seconds




SpO₂ measurement options of Mindray, Masimo SET and Nellcor



Respiration Rate (RR) recorded from SpO₂ finger sensor with Mindray TrueRespTM technology



ARTEMA Sidestream CO₂ with 50ml/min sampling rate

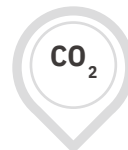

Options of SmarTempTM,™ 3, Genius Exergen TemporalScannerTM

Truly exible work ows

Providing comprehensive documentation at the point of care, the VS 9's workflows unlock the full potential of bedside observations and help consolidate and speed up bedside assessments.

Vivid user interface

- 10.1" high-resolution screen to view more information
- Intuitive gesture operation with capacitive touchscreen

Capture more at the bedside

- Up to 30 customizable manual inputs
- Infographic pain charting tool

Create the right workflow

- Customizable screen layout and parameter settings
- Configurable workflows by department profiles
- Improve caregiver compliance to hospital guidelines

